



Various Functions

■ RX Head (Streak Retinoscope)

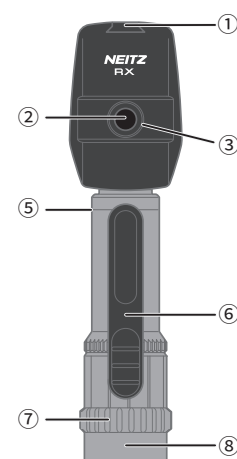


- ① Head Attachment Mount
- ② Viewing Window
- ③ Viewing Window Rim
- ④ Illumination Opening
- ⑤ Reference Line
- ⑥ Streak Adjustment Lever
- ⑦ Streak Rotation Ring
- ⑧ Head Joint

■ RX-SP Head (Spot Retinoscope)

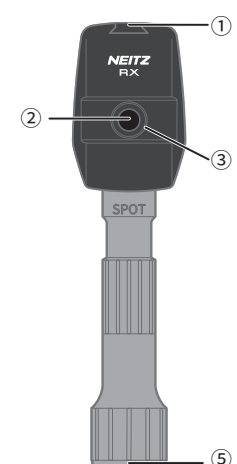


- ① Head Attachment Mount
- ② Viewing Window
- ③ Viewing Window Rim
- ④ Illumination Opening
- ⑤ Head Joint

Charging Stand RC-III

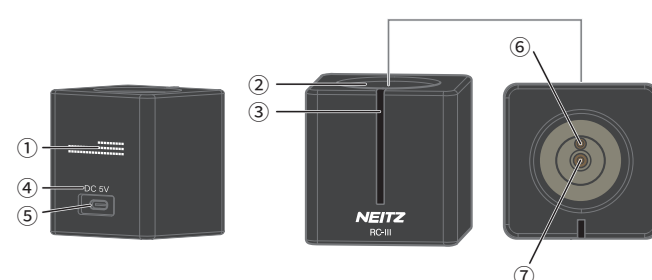

- ① Serial Number and Information Mark
- ② Charging Port
- ③ LED Indicator
- ④ Rated Voltage
- ⑤ USB Socket (USB Type-C)
- ⑥ Negative Contact
- ⑦ Positive Contact

Fixation Card Set FC (Optional)

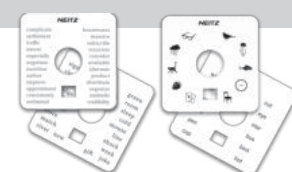

Two Types
of fixation card
(Fixation Card Set RX)



Card Holder RX
(Fixation Card Set RX)



Card Holder RX-SP
(Fixation Card Set RX)

Accessory



Carrying case



Neitz Instruments Co., Ltd.
4F Ichibancho Court, 15-21, Ichibancho,
Chiyoda-ku, Tokyo 102-0082, Japan
Phone: +81-3-3237-0552 Fax: +81-3-3237-0554
URL: <https://www.neitz.co.jp/en/>

Distributor

* As part of our policy of continued product improvement, we reserve the right to alter and/or amend specifications at any time without prior notice.
* The colors in the catalog may be slightly different from those of the actual products .

NEITZ
Made by Neitz for Your Medical Field

RETINOSCOPE
RX



Streak Retinoscope



RX-RC
RX-3
RX-3A

Spot Retinoscope

RX-3SP

Charging Stand

RC-III

The Features Various Functions of Retinoscopes



Smooth Illumination Control

- Smooth stepless light control from the minimum to the maximum.
- Reduces the burden of the patient caused by photophobia.
- Clear decision of "With" reflex, "Against" reflex, and Neutralized reflex.



Powered by Li-ion Battery

- The RX-RC can be stored in the table-top charger while charging.
- *The RX-RC is turned off while charging. No need to operate the power switch.

*For RX-RC



Fixation Card Set FC is Available

- Two types of Fixation Card Set for RX and for RX-SP (Optional)



Basic Functions



Pursuit of Visibility

The anti-reflection filter provides wide and clear field of view. This reduces the sense of fatigue due to a long-time use.



Pursuit of Usability

The streak can be easily set at the best position. The rotation and the width change of the streak are separately controlled. Rotating the streak does not affect the streak width. The lightweight retinoscope head rotates freely.



Sharp Streak Light

The precision cut slit realizes the sharp and clear streak light.



Free Streak Rotation

The streak rotates 360 degrees. You can rotate it from any position and is useful for measurement of the axis of myopia.

Furthermore, the retinoscope is equipped with the Head Attachment Mount for eyeglasses wearers and the feature to attach a presbyopia lens easily.



Fixation Card Set

The Fixation Card Set is available for use in dynamic retinoscopy. Consists of two fixation cards and a card holder. (Optional)

Streak Retinoscope RX-3



- Used with C cell battery handle.
- Package Contents:
 - Retinoscope RX head
 - C cell battery handle
 - Spare bulb
 - +2D Presbyopia lens
 - Head attachment
 - Carrying case
- Bulb: Model No. L-28 (3V, 1.2W)

Streak Retinoscope RX-3A



- Used with AA cell battery handle.
- Package Contents:
 - Retinoscope RX head
 - AA cell battery handle
 - Spare bulb
 - +2D Presbyopia lens
 - Head attachment
 - Carrying case
- Bulb: Model No. L-28 (3V, 1.2W)

Streak Retinoscope RX-RC RENEWED



- Used with a Li-ion battery.
- Charged by the special Charging Stand RC-III (separately sold).
- Package Contents:
 - Retinoscope RX head
 - Rechargeable battery handle
 - Spare bulb
 - +2D Presbyopia lens
 - Head attachment
 - Carrying case
- Bulb: Model No. L-27 (4V, 2.6W)

Spot Retinoscope RX-3SP




- Used with C cell battery handle.
- Delivers parallel light in dot shape.
- Package Contents:
 - Retinoscope RX-SP head
 - C cell battery handle
 - Spare bulb
 - Head attachment
 - Carrying case
- Bulb: Model No. L-05 (3V, 1.5W)

※ **RX-3ASP** used with AA cell battery handle is also available.

Fixation Card Set FC



- Specially designed fixation card set for Neitz retinoscope RX series.

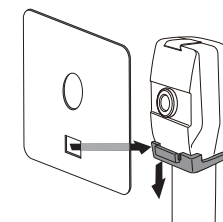

Card Holder

(For RX)

(For RX-SP)

Schematic Eye for Retinoscopy Practice



- The Schematic Eye allows to check the reflex of the "with" and "against" movements by setting one eye at a positive value and the other at a negative value.

- The Schematic Eye shows the difference between the divergent and convergent beam of light emitted from the retinoscope well.

Stand Schematic Eye for Practice

